


# 緊急斷離程序 Emergency Disconnect Procedure


型號 AK 96  
供病人使用  
導管




**1 按血泵的按鈕**  
血泵將停動，按鈕將開始閃動。

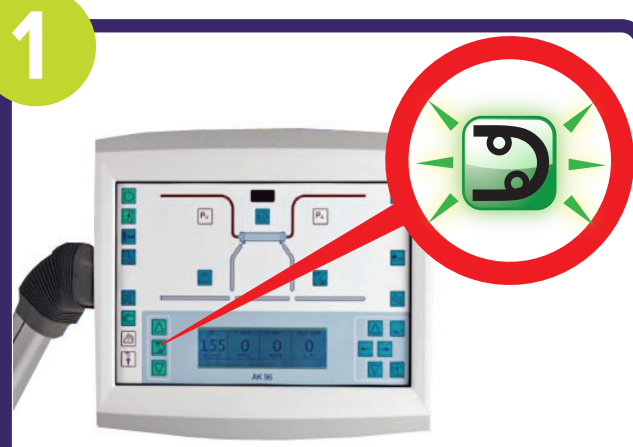

**2 使用在急救袋中  
找到的用品...**

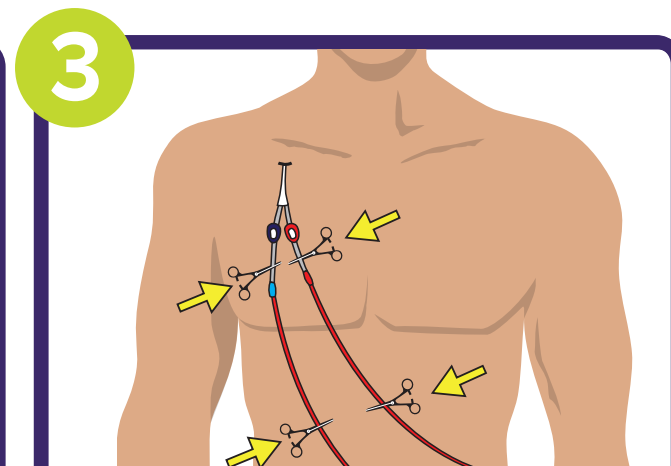

**3 關好導管和血路管上的兩個  
夾子。在連接位下方和上方  
加放安全夾。**

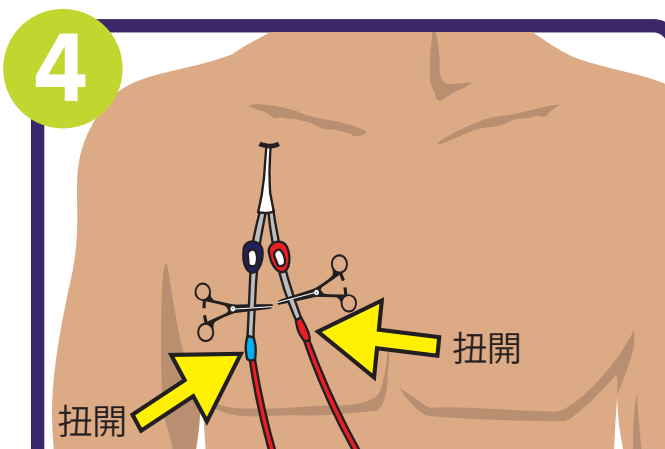

**4 扭開連接位，  
將導管和血路管斷離。**



**5 攜同急救袋  
前往安全的地方**

Produced by **BCRenal**  
Provincial Health Services Authority

**fraserhealth**  
Better health. Best in health care.

**Interior Health**

**northern health**  
the northern way of caring

**Providence  
Health Care**

**Vancouver  
CoastalHealth**  
Promoting wellness. Ensuring care.

**island health**

Provincial Health  
Services Authority

**BRITISH  
COLUMBIA**