


緊急斷離程序 Emergency Disconnect Procedure

型號 AK 96
供病人使用
人工血管或瘻管

1

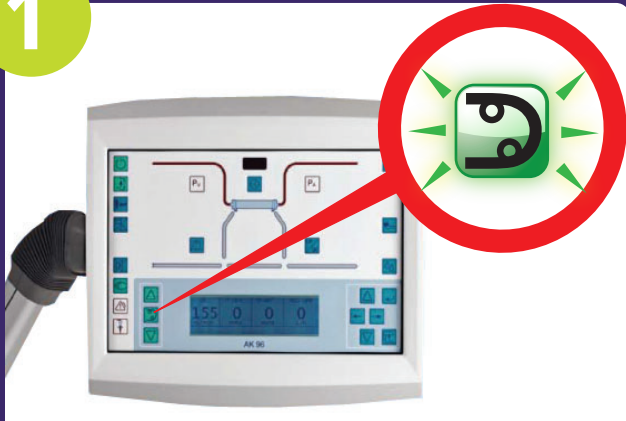


按血泵的按鈕
血泵將停動，按鈕將開始閃動。

2




**使用在急救袋中
找到的用品...**

3

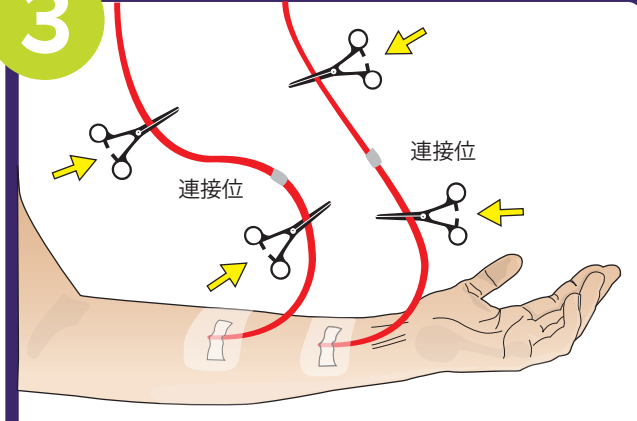

**關好針頭和血路管上的兩個
夾子。在連接位下方和上方
加放安全夾。**

4

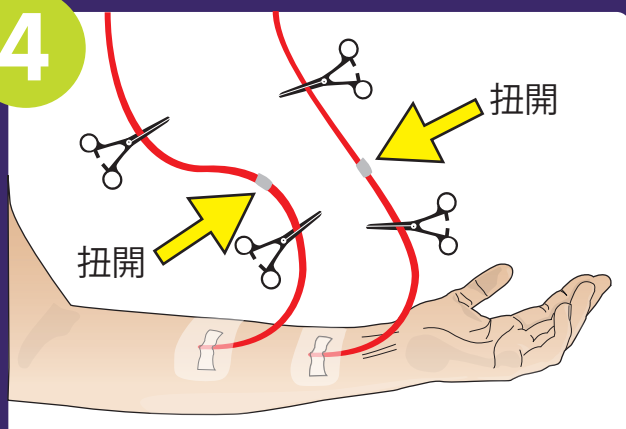

**扭開連接位，
將血路管斷離。**

5



**攜同急救袋
前往安全的地方**

Produced by **BCRenal**
Provincial Health Services Authority

fraserhealth
Better health. Best in health care.

Interior Health

northern health
the northern way of caring

Providence Health Care

Vancouver Coastal Health
Promoting wellness. Ensuring care.

island health

Provincial Health Services Authority

BRITISH COLUMBIA